

TR 4200 Atlas Junior

Bariatric shower trolley



Specifications

Length	2150 mm (85")
Width	740 mm (29")
Height	535-835 mm (21" - 33")
Working Load Limit	270 kg (600 lbs)
Chassis	Powder coated and zinc plated steel
Castors	Brakes on all four castors and straight steering
Control	Hand control for height adjustment and tilting
Mattress, pillow	Flame retardant material
Battery	Rechargeable, 24 V



Registered at Medical Products Agency, Sweden, as a MD Class 1 product, MDD 93/42EEC.

The Atlas Junior Bariatric Shower Trolley - Multi Tasking - Endless Possibilities

The TR 4200 Multi Tasking Shower Trolley is designed for lots of power and loads of space for your most difficult care takers. The trolley is used for shower, transport, as a changing or treatment trolley, the possibilities are endless.

The person can be transferred from bed to shower and back again in a safe and secure manner. The height of the base allows mobile sling lift transfers. Rechargeable battery powered dual actuators raise and lower the stretcher easily and quietly with an easy to use hand control. The hand control operates the power tilt in both directions for improved drainage and tilts to a "Trendelenburg" position when desired. An audible tone sounds when the stretcher is level.

The side rails keep the person safely in place. Side rails swing down for transfer and lock into up position for person security and a secure transport. The trolley has central foot operated castor brakes and straight steering. The stretcher has two positions, flat for transfer and tilt for transport, shower and improved drainage. TR 4200 is delivered with mattress, high capacity rechargeable batteries, hand control, drain plug and hose.



Side rails swing down or lock in horizontal position for transfer and up position for safe transport.



Tilting for "Trendelenburg" position and easy drainage.

TR EQUIPMENT

TR EQUIPMENT AB PO Box 116, S-573 22, TRANÅS, SWEDEN
 TR GROUP INC. 903 WEDEL LANE, GLENVIEW, ILLINOIS 60025, USA
 TR EQUIPMENT UK LTD. INNOVATION CENTRE, GALLOWS HILL, WARWICK CV34 6UW ENGLAND
 TRANSPORHAMED GMBH GARTENSTRASSE 9, 42799 LEICHLINGEN, GERMANY

Phone: +46 140-38 50 60 Fax: +46 140-38 50 61
 Phone: +1 800-752-6900 Fax: +1 847-724-8566
 Phone: +44 844 335 8386 Fax: +44 844 335 8383
 Phone: +49 2175-16 94 22-0 Fax: +49 2175-16 94 22-19

Hygiene Equipment Designed for People

www.trequipment.com



TR 4200 Atlas Junior

Bariatric shower trolley



Features

The stretcher is easily and quietly raised by two powerful actuators. The hand control operates the power tilt in both directions and raises and lowers the trolley.



The stretcher can easily be tilted to the “Trendelenburg” position, to alleviate a person in shock. The height adjustable stretcher enables a comfortable working position for the attendant. An audible tone sounds when stretcher is level. When transporting through narrow doorways the side rails can be folded reducing the trolley width.

The stretcher has rubber bumper wheels on all four corners to protect walls and doorways. Stretcher comes standard with a soft flame retardant mattress with drainage at the foot end and includes a drain hose and pillow.

The trolley is made of durable material, powder coated steel, zinc plated steel that make cleaning simple and efficient. The mattress is easily removable.



Convenient single foot control for straight steering and brakes.

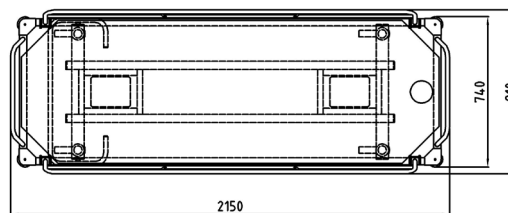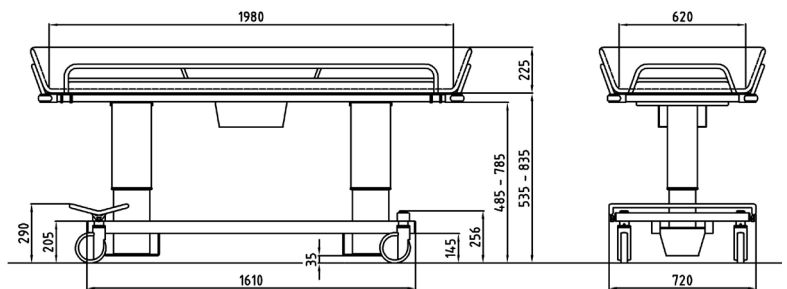

Accessories

- Backrest wedge pillow
- Backrest wedge pillow 15°
- Backrest wedge pillow 30°

TR reserves the right to make design changes and/or improvements without notice

Distributor:

TR EQUIPMENT

TR EQUIPMENT AB PO Box 116, S-573 22, TRANÅS, SWEDEN
 TR GROUP INC. 903 WEDEL LANE, GLENVIEW, ILLINOIS 60025, USA
 TR EQUIPMENT UK LTD. INNOVATION CENTRE, GALLOWS HILL, WARWICK CV34 6UW ENGLAND
 TRANSPORHAMED GMBH GARTENSTRASSE 9, 42799 LEICHLINGEN, GERMANY

Phone: +46 140-38 50 60 Fax: +46 140-38 50 61
 Phone: +1 800-752-6900 Fax: +1 847-724-8566
 Phone: +44 844 335 8386 Fax: +44 844 335 8383
 Phone: +49 2175-16 94 22-0 Fax: +49 2175-16 94 22-19

Hygiene Equipment Designed for People

www.trequipment.com

